




ParaHIT Total+

Rapid Test for *P. falciparum* and Pan Malarial Species

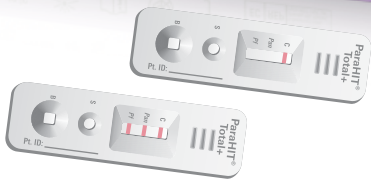

Rapid visual
immunochromatographic test for
detection of **Plasmodial** antigens in whole blood



ParaHIT Total+

An antigen detection assay that detects *P. falciparum* specific Histidine Rich Protein-II (Pf HRP-II) and Pan malarial species specific Plasmodium Lactate Dehydrogenase (Pan pLDH)

Highly Reliable & Accurate Results



Sensitivity (%)	Confidence Interval (% CI)
100	99.01 - 100

Data on Arkray files



Specificity (%)	Confidence Interval (% CI)
99.5	99.34 - 99.95

Data on Arkray files

Convenient and Time Saving




5 µL sample volume, comply with the global product requirement.



Non-instrumental, rapid test procedure of 20 minutes. Suitable for field use.

Flexibility

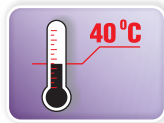

Product is stable up to 40°C, suitable to use in tropical countries.



No sample pre-dilution. Useful for point of care testing & epidemiological studies.

Ordering Information

Product	Sample Transfer Device	Test	25 Test	50 Test
ParaHIT Total+	Plastic Pipette	Device	55IC206-25	55IC206-50

THE MANUFACTURING SITE'S QMS IS
CERTIFIED FOR ISO 13485 & ISO 9001

Arkray Healthcare Pvt. Ltd.

Plot No. 336, 338, 340, Road No. 3, G.I.D.C., Sachin - 394 230, (Surat), INDIA.

Telephone: 0261-239 7712 to 15, Facsimile: 0261-239 7719, Website: www.arkray.co.in

Specifications may change without prior notice